


Orthopaedic Suction Cannula

Product Data Sheet

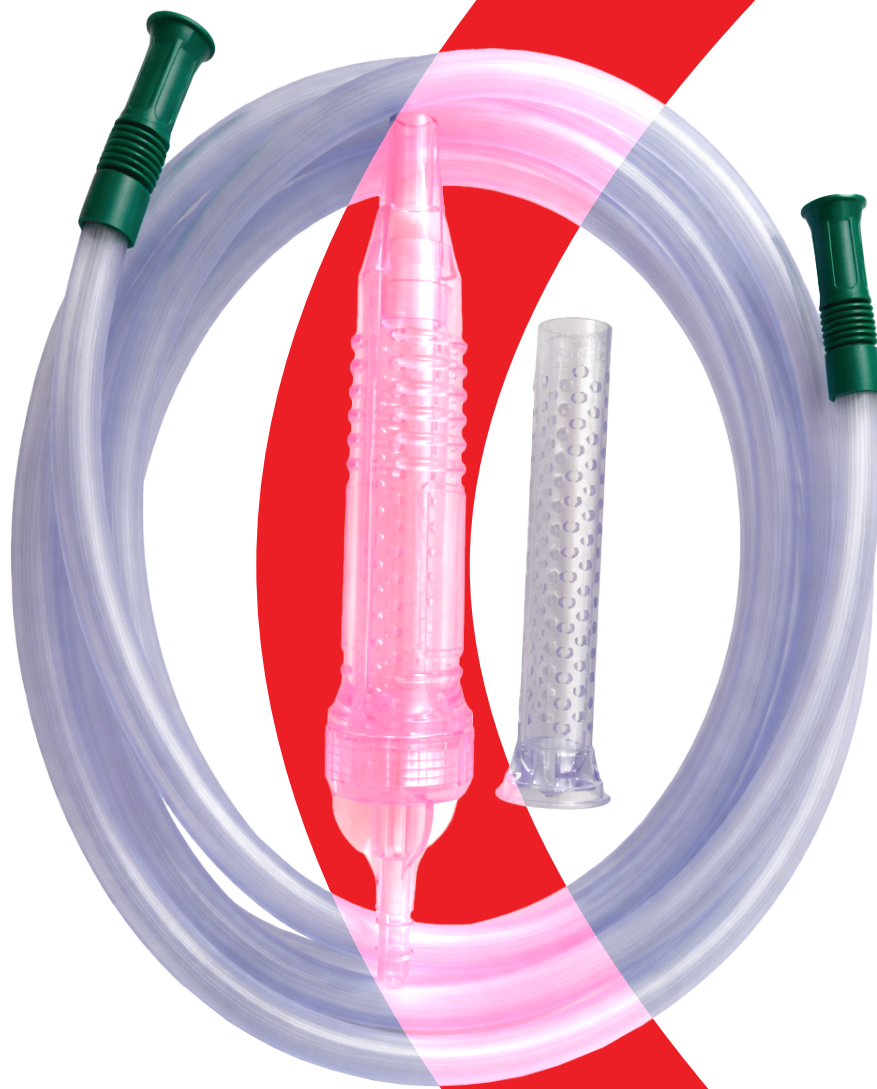

ADVANCEDMEDTECH

Engineered For Tomorrow, For Better Healthcare Today

Orthopaedic Suction Set

Precision. Control. Reliability in the Operating Theatre.

The Orthopaedic Suction Set is a sterile, disposable device. It consists of cannula kits and 3m tubing with connectors, made from medical-grade raw materials. The cannula can be easily disconnected from the suction tube and replaced if necessary. The built-in filter retains bone crumbs inside the cannula while allowing free fluid passage. The filter can be replaced when required. Replaceable tips facilitate adaptation of the cannula for specific procedural conditions.

Intended Use

Orthopaedic suction cannula and Orthopaedic Suction Set are used to draw liquids from the patient during orthopaedic surgery procedures where bone crumbs are expected in the surgical field.

Key Features and Benefits

Efficient Fluid Removal

Maintains a clear surgical field during procedures

Bone Fragment Protection

Prevents blockage caused by bone debris

Rapid Filter Exchange

Enables quick changes without disrupting workflow

Ergonomic Hand Design

Enhances comfort and control for prolonged use

Extended Reach

Reversible suction tip provides extra length when required



Product Specifications

Reference: AMT001

NAPPI Code:

Duration of Use

Intended for continuous use for less than 60 minutes (transient).

Materials

Cannula (handle, filter, chamber): PC

Tube and Connector: PVC

Cannula Tip: PVC

Supplied With

Tube Length: 3 m
Tube Diameter: 7 mm (ID)
Cannula Tip Diameter: 6 mm (ID)
Straight Cannula Tip Length: 25 cm
Curved Cannula Tip Length: 15 cm and 25 cm
xFilter Hole Diameter: 2.5 mm

Clinical Advantages

Suitable for removal of bone fragments during orthopaedic surgery
Eliminates blockages, ensuring uninterrupted suction
Designed for reliability in demanding surgical environments

Packaging and Compliance

Packaging: 6 units per outer carton
Quality System: ISO 13485
EU Classification: Class IIa
Sterilisation Method: Ethylene Oxide

Shelf Life and Storage

Shelf Life: 5 years
Storage Conditions: Store in a clean, dry place away from direct sunlight



Manufactured for
ADVANCEDMEDTECH

Address: Unit 12, Paramount Park,
Mount Edgecombe, Durban, 4003